

***In vitro* Assay for Dendritic Spine Retraction of Hippocampal Neurons with Sparse Labeling**

Ryohei Iwata^{1, 2, 3, *} and Takuji Iwasato^{1, 2}

¹Division of Neurogenetics, National Institute of Genetics, Mishima, Japan; ²Department of Genetics, SOKENDAI, Mishima, Japan; ³ULB Neuroscience Institute (UNI), Institut de Recherches en Biologie Humaine et Moleculaire (IRIBHM), Universite Libre de Bruxelles (ULB), Brussels, Belgium

*For correspondence: ryohei.iwata@ulb.ac.be

[Abstract] Dendritic spines are the post-synaptic structures that play a central role in excitatory synaptic transmission. Developmental spinogenesis relies on a variety of stimuli such as those derived from cell-cell communication and their downstream signaling. Here, we describe an *in vitro* assay of dendritic spine retraction using hippocampal slice culture, in which individual neurons are sparsely and brightly labeled by the Supernova method, for the study of molecular mechanisms of spine development.

Materials and Reagents

A. Materials

1. Electrode (5 mm Φ platinum disk) (Nepa Gene, model: CUY650P5)
2. Surgical needle (ELP, model: CR13-50)
3. Glass capillary (Warner Instruments, model: GC150TF-10)
4. NuncTM culture plate (6-well) (Thermo Fisher Scientific, Thermo ScientificTM, catalog number: 140675)
5. Millicell culture plate inserts (EDM Millipore, catalog number: PICM03050)
6. Confetti (LCR membrane filter) (EDM Millipore, catalog number: FHLC01300)
7. Petri dish (60 mm) (Sigma-Aldrich, catalog number: Z721034)
Note: This product has been discontinued.
8. Microcentrifuge tube (Sigma-Aldrich, catalog number: Z666505)
9. Cover glass (MATSUNAMI GLASS)
10. Microscope glass slide (MATSUNAMI GLASS)
11. Aspirator tube assembly (Drummond Scientific, catalog number: 2-000-000)

B. Sparse and bright labeling of hippocampal neurons using the *in utero* electroporation-based Supernova system (Mizuno *et al.*, 2014)

1. Pregnant mouse (E13.5-E15.5)
2. Maxi-prep kit
3. Trypan blue solution (0.4%) (Sigma-Aldrich, catalog number: T8154)
4. Supernova vector DNA solution (pK036.TRE-FIpe-WPRE-pA and pK037.CAG-FRT-STOP-FRT-TurboRFP-ires-tTA-WPRE-pA)

5. Somnopentyl (Pentobarbital sodium)
6. 70% ethanol
7. Saline

C. Organotypic hippocampal slice culture

1. Mouse pup (P4-P5)
2. MEM with gluramax-1 (Thermo Fisher Scientific, catalog number: 41090-028)
3. EBSS (Thermo Fisher Scientific, catalog number: 14155-048)
4. D-glucose (NACALAI TESQUE, catalog number: 16805-35)
5. Penicillin-streptomycin (Thermo Fisher Scientific, Gibco™, catalog number: 15140-122)
6. Nystatin (Thermo Fisher Scientific, Gibco™, catalog number: 15340029)
7. HEPES (1 M) (Thermo Fisher Scientific, Gibco™, catalog number: 15630-080)
8. Horse serum (Sigma-Aldrich, catalog number: H1270-500 ml)
9. Culture medium (See Recipes)
10. Slicing buffer(See Recipes)

D. *In vitro* retraction assay

1. EphrinA3-Fc (R&D Systems, catalog number: BT359)
2. Human Fc fragment (Jackson ImmunoResearch, catalog number: 009-000-008)
3. Goat anti-human IgG, Fc fragment specific (Jackson ImmunoResearch, catalog number: 109-001-008)
4. Paraformaldehyde (NACALAI TESQUE, catalog number: 02890-45)
5. PBS (-)
6. Vectashield mounting media with DAPI (Vector laboratories, catalog number: H-1200)

Equipment

A. Sparse and bright labeling of hippocampal neurons using the *in utero* electroporation-based Supernova system.

1. Electroporator (Nepa Gene, model: Nepa21)
2. Water bath
3. Puller (NARISHIGE Group, model: PC-10)
4. Heating table (Leica Biosystems, model: HI1220)
5. Surgical tools (scissor, 12 cm) (Fine Science Tools, catalog number: 14001-12)
6. Surgical tools (scissor, 10.5 cm) (Fine Science Tools, catalog number: 14088-10)
7. Surgical tools (forceps, 12 cm) (Fine Science Tools, catalog number: 11000-12)
8. Surgical tools (forceps, 10 cm) (Fine Science Tools, catalog number: 11050-10)
9. Surgical tools (ring forceps, 8.5 cm) (Fine Science Tools, catalog number: 11101-09)

- B. Organotypic hippocampal slice culture
 - 1. LinearSlicer (Dosaka, model: PRO7)
 - 2. 37 °C, 5% CO₂ incubator
 - 3. Pipette
- C. *In vitro* retraction assay
 - 1. 25 °C incubator
- D. Imaging and analysis
 - 1. Confocal microscope (Leica Biosystems, model: TSC-SP5)
 - 2. Imaris filament tracer (Bitplane)

Software

- 1. Imaris

Procedure

- A. Sparse and bright labeling of hippocampal neurons using the *in utero* electroporation-based Supernova system

Note: Using the Supernova system that enables sparse and bright cell labeling with little background is highly recommended for analyses of fine structures such as spines of individual neurons. The Supernova system is described elsewhere (Mizuno et al., 2014). Briefly, this system uses two vectors, TRE-Flpe-WPRE-pA (pK036) and CAG-FRT-STOP-FRT-TurboRFP-ires-tTA-WPRE-pA (pK037). In a small population of neurons that carry both vectors, TRE leakage drives weak but over-threshold Flpe expression from the first vector, which is followed by removal of the FRT-STOP-FRT cassette in a few copies of the second vector and weak expression of tTA from these copies. Then, only in these sparse neurons, tTA binds with TRE and induces strong Flpe expression, which results in removal of the FRT-STOP-FRT cassette from many copies of the second vector, and finally results in extremely strong RFP expression via the positive feedback. Flpe-based Supernova vectors [TRE-Flpe-WPRE-pA (pK036) and CAG-FRT-STOP-FRT-TurboRFP-ires-tTA-WPRE-pA (pK037)] are available from T.I. (tiwasato@nig.ac.jp).

- 1. To prepare micropipettes, pull glass capillary using a puller under the following conditions: two-step; heater, 65-75; weight, 1-3. Cut off the tip with forceps for appropriate diameter (about 20 µm).
- 2. To prepare DNA solution for the Supernova labeling, purify the plasmid set, TRE-Flpe-WPRE-pA (pK036) and CAG-FRT-STOP-FRT-TurboRFP-ires-tTA-WPRE-pA (pK037), using endotoxin free Maxi-prep kit, and mix them to a final concentration of 1 µg/µl of pK036 and 50 ng/µl of pK037.

3. Aliquot 40 μ l of the DNA solution, and add 4 μ l of Trypan blue solution (0.4%).
4. Anesthetize a pregnant mouse (E13.5-E15.5) with pentobarbital sodium through intraperitoneal injection (50 mg/kg in saline).
5. Place the anesthetized mouse on a working plate, and wash the abdomen with 70% ethanol.
6. Make an incision at the abdominal midline using scissors, and take out the uterus gently onto the saline-moistened gauze. During the surgery, keep the uterus moist with pre-warmed saline (Figure 1).

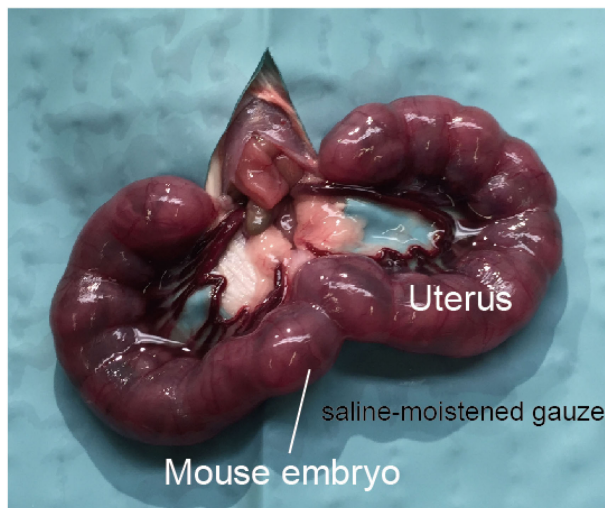


Figure 1. An anesthetized pregnant mouse with embryos (E14.5)

7. Inject 1-2 μ l of the DNA solution into the target ventricle of the embryos using the micropipette with an aspirator tube.
8. Hold the embryos with electrodes (Figure 2: negative electrode on the injected side of the embryo head), and apply electric pulses [40-50 V, P on 50 msec, P off 950 msec (1 Hz), 5 times].

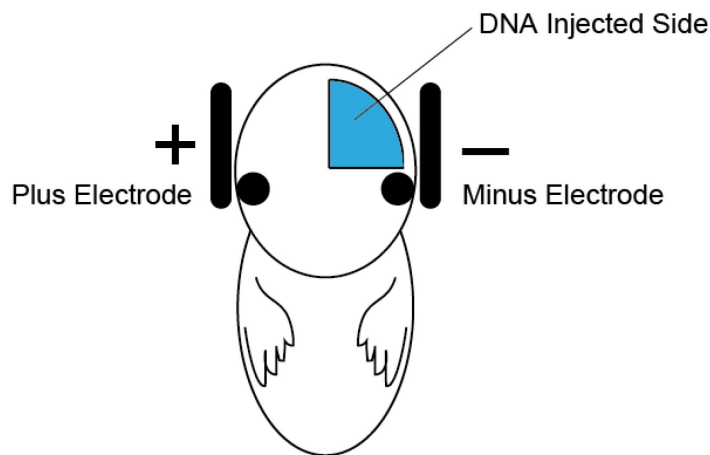


Figure 2. Schematic illustration of *in utero* electroporation. DNA injection side (blue) and the relative position of the electrodes for the hippocampal electroporation.

9. Place the uterus gently into the abdominal cavity.
10. Fill the cavity with pre-warmed saline.
11. Suture the surgical incision using surgical needle and thread.
12. Place the mouse on a heating table at 37 °C until recovery from anesthesia.

B. Organotypic hippocampal slice culture

1. Pour 1 ml of culture medium into each well (6-well plate).
2. Place a Millicell culture plate insert into each well.
3. Cut confetti disc using scissors (quarter), and put a piece on each insert.
4. Place the plate into a 5% CO₂ incubator at 37 °C.
5. Pour 10 ml of slicing buffer into a Falcon tube (50 ml), and keep it on ice.
6. Pour 5 ml of slicing buffer into a Petri dish (60 mm) on ice.
7. Euthanize pup (P4-P5) by decapitation using scissors.
8. Place the brain into a Petri dish containing slicing buffer on ice.
9. Cut the brain half along the midline after removing the cerebellum.
10. Glue the brain tissue on the stage of LinearSlicer buffer tray.
11. Fill the buffer tray with the slicing buffer (4 °C).
12. Cut a 300 µm hippocampal slices (sagittal sections) using LinearSlicer. Leave the entorhinal cortex attached to the hippocampus (Figure 3).

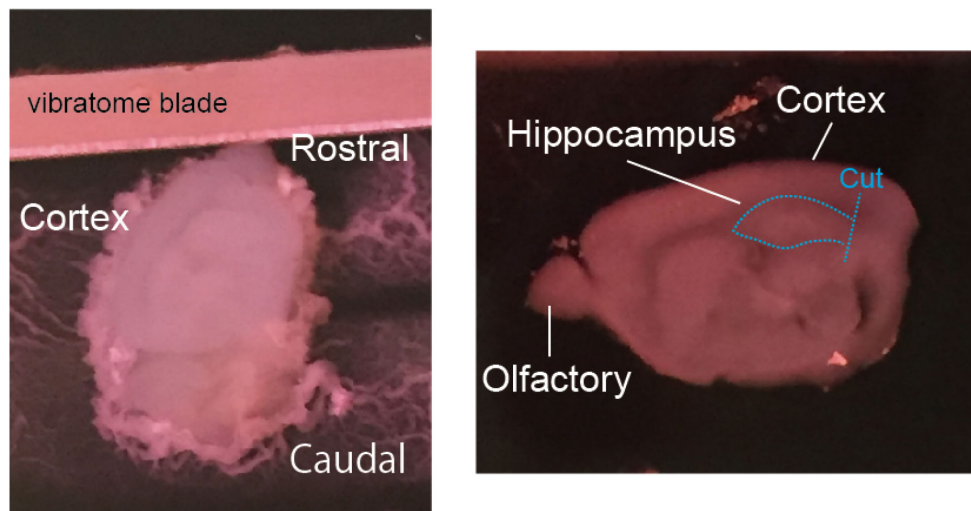


Figure 3. Dissection of the hippocampus from the brain. Left: A P5 mouse brain was sagittally sliced by a vibratome. Right: A 300 µm-thick slice and cut position were shown.

13. Collect the slices in a 50 ml Falcon tube containing slicing buffer (4 °C).
14. Take a slice from the tube using a pipette, and put it on the center of the confetti on insert.
15. Remove excess buffer from the insert, and place the plate in a 5% CO₂ incubator (37 °C).
16. Change culture medium every 2 days until retraction assay.

C. *In vitro* retraction assay

1. Mix 10 µg of ephrinA3-Fc/Fc and 45 µg of anti-human Fc antibody in a microcentrifuge tube.
2. Place the tube into a 25 °C incubator.
3. After 1 h, place the pre-clustered solution into organotypic hippocampal slice culture medium.
The final concentration for ephrinA3-Fc/Fc is 10 µg/ml.
4. Place the plate into a 5% CO₂ incubator at 37 °C.
5. After 16 h, fix the hippocampal slices on the confetti discs with 4% PFA in PBS for 30 min at room temperature.
6. Mount hippocampal slices onto a glass slide with Vectashield mounting media, and cover with a cover glass.

D. Imaging and analysis

1. Image dendritic spines in primary/secondary dendrites of hippocampal pyramidal neurons in CA1 region using a confocal microscope. Sequential z-images consisted of optical section (1,024 x 1,024 pixels) with 0.1 µm intervals using 63x oil immersion objective (numerical aperture, 1.3) with 9x digital zoom. Figure 4 shows z-stacked images of typical dendritic spines.
2. Reconstruct z-images, and measure spine density, spine length, volume by the Imaris Filament Tracer.

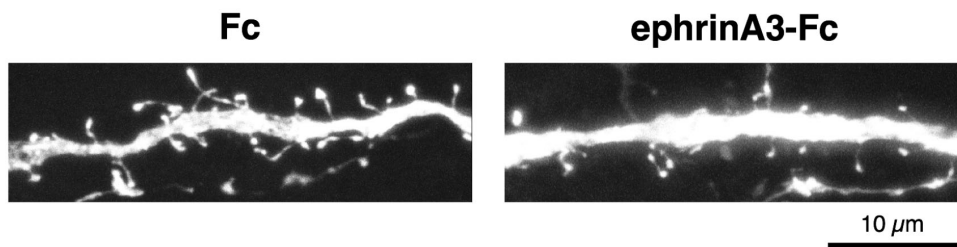


Figure 4. Example of results of dendritic spine retraction assay. Segments of dendrites from hippocampal slices, after treatment with ephrinA3-Fc (right) or Fc (left). Scale bar: 10 μ m.

Recipes

1. Culture medium
250 ml MEM with glutamax-1
125 ml horse serum
120 ml EBSS with 3% D-glucose
5 ml penicillin-streptomycin
0.3 ml nystatin
2. Slicing buffer
487.5 ml EBSS
12.5 ml 1 M HEPES

Acknowledgments

This protocol was adopted from Iwata *et al.* (2015). This work was supported by the JSPS KAKENHI (15H04263, 16K14559) and MEXT KAKENHI (15H01454).

References

1. Iwata, R., Matsukawa, H., Yasuda, K., Mizuno, H., Itohara, S. and Iwasato, T. (2015). [Developmental RacGAP \$\alpha\$ 2-chimaerin signaling is a determinant of the morphological features of dendritic spines in adulthood.](#) *J Neurosci* 35(40): 13728-13744.
2. Mizuno, H., Luo, W., Tarusawa, E., Saito, Y. M., Sato, T., Yoshimura, Y., Itohara, S. and Iwasato, T. (2014). [NMDAR-regulated dynamics of layer 4 neuronal dendrites during thalamocortical reorganization in neonates.](#) *Neuron* 82(2): 365-379.