

Axonal Conduction Velocity Measurement

Margaret Louise DeMaegd, Carola Städele* and Wolfgang Stein*

School of Biological Sciences, Illinois State University, Normal, IL, USA

*For correspondence: carola@neurobiologie.de; wstein@neurobiologie.de

[Abstract] Action potential conduction velocity is the speed at which an action potential (AP) propagates along an axon. Measuring AP conduction velocity is instrumental in determining neuron health, function, and computational capability, as well as in determining short-term dynamics of neuronal communication and AP initiation (Ballo and Bucher, 2009; Bullock, 1951; Meeks and Mennerick, 2007; Rosenthal and Bezanilla, 2000; Städele and Stein, 2016; Swadlow and Waxman, 1976). Conduction velocity can be measured using extracellular recordings along the nerve through which the axon projects. Depending on the number of axons in the nerve, AP velocities of individual or many axons can be detected.

This protocol outlines how to measure AP conduction velocity of (A) stimulated APs and (B) spontaneously generated APs by using two spatially distant extracellular electrodes. Although an invertebrate nervous system is used here, the principles of this technique are universal and can be easily adjusted to other nervous system preparations (including vertebrates).

Keywords: Nervous system, Action potential, Propagation, Extracellular stimulation, Ectopic, Nerve, Spike timing, Axon

[Background] Long-distance communication in the nervous system is mediated by APs that travel along axons. The ionic currents that flow across the axon membrane when an AP is generated (Hodgkin and Huxley, 1952) can be detected even outside of the neuron, using extracellular recording electrodes. AP conduction velocities in different neurons are quite variable and range from 200 meters per second (447 miles per hour) to less than 0.1 meters per second (0.2 miles per hour) (Kress *et al.*, 2008; Kusano, 1966). In order to understand why there are differences in conduction velocity, the passive (membrane) properties of the axon need to be taken into account. Some axons propagate information more rapidly than others because of differences in properties that affect the time constant (e.g., resistance and capacitance) and the length constant (e.g., axon diameter, membrane permeability, and degree of myelination). Especially in unmyelinated axons, conduction velocity largely depends on the axon diameter, which in turn is also correlated with the amplitude of the extracellular AP (Stein and Pearson, 1971). Consequently, determining AP conduction velocity provides more than just information about signal movement and timing. It can also be used to characterize changes in intrinsic axon properties.

Materials and Reagents

Note: The materials and equipment listed refer to the equipment used in Städele and Stein (2016). To reduce costs, comparable materials, equipment and software may be used that serve the same functions. For the general public or a teaching classroom we suggest utilizing equipment from Backyard Brains (<http://backyardbrains.com>).

1. Petri dish (100 x 15 mm, Fisher Scientific, catalog number: FB0875713) lined with silicon elastomer (e.g., Sylgard 184, Sigma-Aldrich, catalog number: 761036; or Elastosil RT 601, Wacker Chemie, catalog number: 60063613)
2. Minutien pins (Fine Science Tools, catalog number: 26002-10)
3. Modeling clay (craft store)
4. Syringe, filled with petroleum jelly (100% pure, pharmacy)
5. Recording/stimulation electrodes

Note: For details how to prepare the petroleum jelly filled syringe or recording/stimulating electrodes, see our companion protocol 'Extracellular axon stimulation' by Städele, C., DeMaegd, M. and Stein, W.

6. Dissected nervous system

*Note: We are using adult Jonah crabs (*Cancer borealis*), purchased from The Fresh Lobster Company (Gloucester, MA).*

7. Physiological saline (see Recipes)

The recipe for *C. borealis* saline can be found in Table 1

- a. Sodium chloride, NaCl (Sigma-Aldrich, catalog number: S9625)
- b. Magnesium chloride hexahydrate, $\text{MgCl}_2 \cdot 6\text{H}_2\text{O}$ (Sigma-Aldrich, catalog number: M9272)
- c. Calcium chloride dihydrate, $\text{CaCl}_2 \cdot 2\text{H}_2\text{O}$ (Sigma-Aldrich, catalog number: C7902)
- d. Potassium chloride, KCl (Sigma-Aldrich, catalog number: P9541)
- e. Trizma base (Sigma-Aldrich, catalog number: T1503)
- f. Maleic acid (Sigma-Aldrich, catalog number: M0375)

Equipment

1. Ruler or micrometer scale
2. Stereomicroscope (e.g., Leica Microsystems, model: MS5)
3. Stimulator (A.M.P.I, model: Master 8 Pulse Stimulator)
Low-cost alternative: Pulse Pal V2 (Sanworks, catalog number: 1102)
4. Amplifier (A-M Systems, model: Differential AC Amplifier 1700, catalog number: 690000)
Low-cost alternative: Spikerbox (Backyard Brains, model: Neuron SpikerBox)
5. Data acquisition board (Cambridge Electronic Design Limited, model: Power 1401-3A)

Low-cost alternative: by using the BYB Spike Recorder, data can be digitized by using the microphone jack and soundcard on a computer/laptop. A second low-cost alternative is Spikehound (<http://spikehound.sourceforge.net>), which also allows recording through the computer soundcard.

6. Camera (e.g., AmScope, model: 3MP USB2.0 Microscope Digital Camera, catalog number MU300)

Software

1. ImageJ (National Institute of Health) or a comparable software
2. Recording software (Spike2 version 7.12, Cambridge Electronic Design Limited)

Low-cost alternative: BYB Spike Recorder (freeware, available on <https://backyardbrains.com/products/spikerecorder>) or Spike hound (<http://spikehound.sourceforge.net>)

Procedure

Figure 1 illustrates the experimental setup and electrode placement. In order to determine AP conduction velocity, it is important to keep in mind that there are differences in the procedure for measuring stimulated APs (Figure 1A) and spontaneously generated APs (Figure 1B). Especially in the latter case, it is crucial that APs are not initiated between the two recording sites, as this will result in wrong velocity measurements. APs can spread bi-directionally from their site of initiation, and if APs are initiated between the two recording sites, the measured difference in AP arrival times between recording sites is not representative of APs travelling the distance between the two (Figure 1C, and also step B2). Thus, in the following, detailed procedures are given for measuring conduction velocity of:

- A. Stimulated APs by means of extracellular stimulation
- B. Spontaneously initiated APs

In both cases, AP conduction velocity is measured as a function of the time the AP needs to travel from one site on an axon/nerve to another site (= 'latency'). Conduction velocity can be calculated using:


Equation 1:
$$v = \frac{\Delta x}{\Delta t}$$

With v as the conduction velocity,


Δx as the spatial distance between recording sites,

Δt as the time difference between recording sites (latency).

A. Stimulated APs

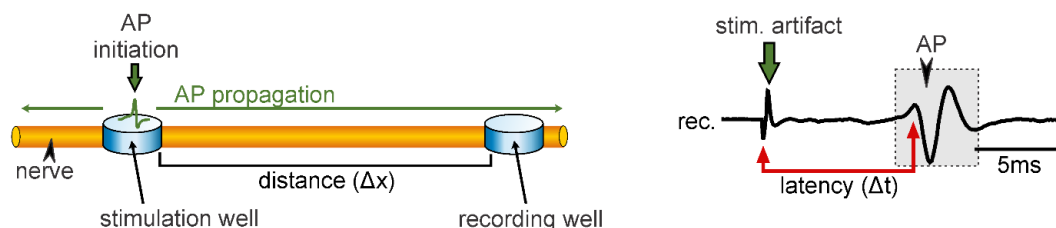

B. Spontaneously generated APs

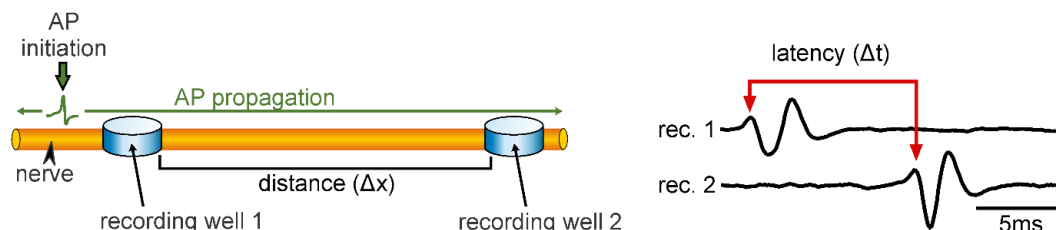

C. Spontaneously generated APs (between recording wells)

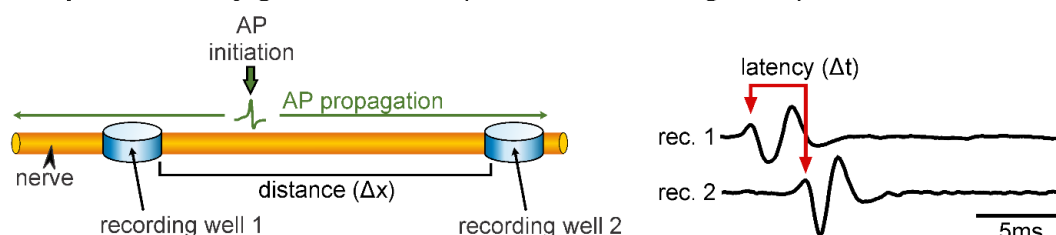


Figure 1. Representation of the electrode placement for measuring conduction velocity of stimulated APs (A) and spontaneously generated APs (B & C). A. Left: APs are initiated within the stimulation well and propagate along the axon (green arrows) towards the recording well. Note that APs propagate bi-directionally but are only detected at the recording site. Δx represents the distance between the stimulation and recording wells. Right: Extracellular recording of a nerve showing the latency (Δt) between the onset of the stimulation (stim. artifact) and the AP (gray box). B. Left: Correct electrode placement for measuring conduction velocity of spontaneously generated APs. Recording wells need to be placed in such a way that the AP is passing the wells without being initiated in-between the wells. Right: The latency in AP occurrence at the two recording sites reflects the time the AP needs to propagate from recording 1 to recording 2. C. Left: Incorrect electrode placement for measuring conduction velocity of spontaneously generated APs. Right: The latency between the two recording sites is small and it only reflects that the AP was initiated spatially closer to recording 1 than to recording 2.

A. Measuring conduction velocity of stimulated APs

The advantage of measuring conduction velocity of stimulated APs is that AP initiation is controlled by the experimenter. Thus, different AP frequencies can be tested, and APs can reliably be elicited

in many different conditions (e.g., at different temperatures or in distinct neuromodulatory conditions). However, because the largest axons respond to the lowest stimulation amplitudes, measuring the velocities of small diameter axons may not be feasible with extracellular stimulation.

1. Dissect the nerve of interest.

*Note: The dissected nerve should be long enough that the two wells constructed in step 2 are separate, and the entire length of the nerve contains the axon of interest. For a detailed protocol how to dissect the stomatogastric nervous system of *C. borealis* please review Gutierrez and Grashow (2009).*

2. Establish the extracellular recording and stimulation as described in companion protocol 'Extracellular axon stimulation' by Städele, C., DeMaegd, M. and Stein, W. For electrode placement see Figure 1A.
3. After determining the AP stimulation threshold, continuously stimulate the axon with a low frequency stimulation (1 to 2 Hz).
4. Identify the AP on the extracellular recording (Figure 1A, right).

Note: The stimulated AP will travel from the stimulation electrode to the recording site. From one stimulus pulse to the next, AP amplitude and shape should not change. Similarly, the latency between the stimulation onset and the AP occurrence should remain constant. If there is no consistency of AP amplitude or latency, then either no AP was elicited or the elicited AP does not pass the recording site, indicating that the axon might be damaged or might project through a different nerve.

5. Measure conduction velocity as described in Data analysis.
6. Variation of the experiment: In many axons, AP conduction velocity is frequency-dependent (Miller and Rinzel, 1981; Weidner *et al.*, 2002) and/or history-dependent (Ballo *et al.*, 2012; Kiernan *et al.*, 1997). To test whether this is the case for the axon of interest, increase the stimulation frequency and observe whether conduction velocity changes as a function of frequency and/or number of APs. Instead of stimulating at a constant frequency, more biological relevant pulse sequences may be applied, such as burst-like changes in frequency (Ballo and Bucher, 2009).

B. Measuring conduction velocity of spontaneously generated APs

Conduction velocity measurements of spontaneously initiated APs does not require extracellular nerve stimulation. Many axons spontaneously generate APs, and accordingly, most nerve recordings show spontaneously occurring APs. By using spontaneously generated APs, the conduction velocities of small diameter axons can be assessed. However, the occurrence and frequency of the APs cannot be controlled by the experimenter. We are using the anterior gastric receptor neuron (AGR; Smarandache and Stein, 2007) of the stomatogastric nervous system of the Jonah crab (*Cancer borealis*) as an example. AGR spontaneously generates tonic spike activity (2 to 5 Hz) in its axon trunk that can be measured with extracellular electrodes (Daur *et al.*, 2009; Städele and Stein, 2016).

1. Dissect the nerve of interest. Then, instead of one stimulation and one recording well, both will be used as recording wells (Figure 1B).
2. Connect the cable of the stimulation electrode to a second channel of the extracellular amplifier.
Note: If the neuron of interest is bipolar or able to generate ectopic APs (e.g., Daur et al., 2009), make sure that the site where APs are initiated is not in between the two recording wells since the measured latencies will not reflect travel times between recording sites (Figure 1C). If necessary, this can be determined by adding a third recording well. The order at which the APs arrive at the recordings wells will be indicative of the initiation site. If either of the two outside wells records the AP first, the correct velocity can be measured using the difference in arrival times between the other recording wells. If the center well records the AP first, use longest latency between center and outside wells.
3. Identify the AP of interest on the extracellular recordings.
Note: The AP of interest should appear with a constant latency between both recording sites. Its shape and amplitude may differ between recording sites, but should be consistent for each site. If the AP of interest is not the largest unit on the extracellular recording, or difficult to identify for other reasons, see Notes.
4. Measure conduction velocity as described in Data analysis.

Data analysis

1. Measure the latency between AP occurrence at the two recording sites.
 - a. For stimulated APs, measure the time difference between stimulus and recorded AP, as illustrated in Figure 1A.
Note: Stimulation onset can be determined using the stimulus artifact, which spreads almost instantaneously through the Petri dish. In the recording software, measure the time between the peak of the stimulus artifact and the peak of the recorded action potential.
 - b. For spontaneously elicited APs, measure the time difference between the peaks of the AP at recording 1 and recording 2, as illustrated in Figure 1B.
2. Measure the distance between recording sites.
Use a ruler or micrometer scale to measure the length of the nerve between the stimulation and recording wells. The shortest distance between both wells should be measured (well edge to well edge). There are different methods that can be used to measure distance. The simplest and least accurate method is using a ruler or micrometer. A more accurate method is the following:
 - a. Put a size standard underneath the Petri dish (e.g., micrometer scale).
 - b. Take a photo that shows the distance between the stimulation and recording well plus the size standard.

- c. Use ImageJ (National Institute of Health) or a comparable software to measure the nerve distance between the two wells. Use the size standard to calibrate the distance measurement in the software.
3. Calculate the conduction velocity using Equation 1.
Note: For more accurate results, calculate, then average the conduction velocities of at least 10 APs.
4. Use a paired *t*-test when comparing conduction velocities in two different conditions (e.g., before and after neuromodulator application) from the same preparation with a significance cut off value of $P < 0.05$.
5. Use a repeated measures analysis of variance (ANOVA) when comparing multiple measurements of conduction velocity in different conditions (e.g., application of multiple neuromodulators, or measurements are multiple temperatures) from the same animal

Notes

1. Using multisweeps and averages to increase measurement accuracy. To measure the latency between stimulus artifact and recording (or between recording 1 and recording 2), it is helpful to not only use single APs. Rather, multiple APs should be taken into account (Video 1). This will not only make the identification of time-correlated APs easier, but also allows to assess variability and changes in conduction velocity that may occur in different conditions. Time-correlated APs are most easily detected using a triggered sweep display (oscilloscope view, or 'multisweep', Video 2). In this view, multiple occurrences of the same AP are superimposed to show APs that occur at a fixed latency to the stimulation or the other recording, respectively. For this, an amplitude threshold is set in the software, through which the AP must rise to start a new sweep ('triggering'). Consequently, only APs large enough to cross the threshold trigger new sweeps, and the latency of arrival of these APs at the recording site becomes obvious. If the threshold is set lower, APs of multiple axons may be detected. Depending on their conduction velocity, they may arrive at distinct, but fixed times at the recording. The APs of other axons will not be time locked to the trigger and hence occur randomly on the display. Using a triggered sweep display is particularly helpful to detect spontaneously generated APs, since it is difficult to determine time correlations between two recording sites using individual APs or without proper analysis.

To determine the exact latency, it can be helpful to average several sweeps of triggered APs (Video 2). In this case, all randomly occurring APs will be small or absent, and only the time-locked APs show a clear peak that can easily be measured.

Video 1. Video clip illustrating time-correlated APs using Spike recorder software (Backyard Brains) for spontaneously generated APs. A similar procedure can be used for stimulated APs.



Video 2. Video clip illustrating spike-triggered multisweeps and averaging using Spike 2 software (CED) for spontaneously generated APs. A similar procedure can be used for stimulated APs.



2. Noise issues. Electrical noise is a common problem associated with the recording of neuronal signals. The major source of electrical noise is pickup, which is caused by radiation from electrical devices and the recording equipment. This radiation produces currents in the electrodes and wires leading to the amplifiers in the recording system. Since most devices are powered by 60 Hz (US)/50 Hz (most other countries) alternating current, pickup appears as a distorted sine wave with a period of 16.7 (US)/20 (others) milliseconds. Noise may be reduced by using a grounded, metal enclosure (Faraday cage) around the preparation and the electrodes. The enclosure separates the source of the radiation from the electrodes. Grounding all recording and stimulation equipment to a single source will reduce noise. This can be done using a banana plug with a crocodile clamp on one side to connect the ground ports of the recording and

stimulation equipment to the Faraday cage. For more information on the principles of Faraday cage function, see Chapman *et al.* (2015) and DeFreitas *et al.* (2012). The person operating the equipment might also be a source of noise, and he or she may need to be grounded as well. Additionally, shielded cables may be used to carry the signals from the electrodes to the amplifier. If the signal-to-noise ratio is still small, the issue might originate from a small signal rather than large noise. In this case, the recording well is likely leaky and should be remade.

Recipes

1. *Cancer borealis* physiological saline recipe (see Table 1)

Table 1. *Cancer borealis* physiological saline recipe

Salt (Molecular Weight)	Concentration [mM]
NaCl (58.44)	440
MgCl ₂ ·6H ₂ O (203.31)	26
CaCl ₂ ·2H ₂ O (147.01)	13
KCl (74.56)	11
Trizma base (121.14)	11.2
Maleic acid (283.3)	5.1

Note: Adjust pH to 7.4-7.6 with Trizma base and maleic acid.

Acknowledgments

This work was supported by grants from the Deutsche Forschungsgemeinschaft (DFG STE 937/9-1), National Science Foundation (NSF IOS 1354932), Illinois State University, and the German Academic Exchange Service.

References

1. Ballo, A. W. and Bucher, D. (2009). [Complex intrinsic membrane properties and dopamine shape spiking activity in a motor axon](#). *J Neurosci* 29(16): 5062-5074.
2. Ballo, A. W., Nadim, F. and Bucher, D. (2012). [Dopamine modulation of Ih improves temporal fidelity of spike propagation in an unmyelinated axon](#). *J Neurosci* 32(15): 5106-5119.
3. Bullock, T. H. (1951). [Conduction and transmission of nerve impulses](#). *Annu Rev Physiol* 13: 261-280.
4. Chapman, S. J., Hewett, D. P., and Trefethen, L. N. (2015). [Mathematics of the Faraday cage](#). *Siam Review* 57(3): 398-417.
5. Daur, N., Nadim, F. and Stein, W. (2009). [Regulation of motor patterns by the central spike-initiation zone of a sensory neuron](#). *Eur J Neurosci* 30(5): 808-822.
6. DeFreitas, J. M., Beck, T. W., and Stock, M. S. (2012). [Comparison of methods for removing electromagnetic noise from electromyographic signals](#). *Physiol Meas* 33(2): 147-159.

7. Gutierrez, G. J. and Grashow, R. G. (2009). [Cancer borealis stomatogastric nervous system dissection](#). *J Vis Exp* (25).
8. Hodgkin, A. L. and Huxley, A. F. (1952). [A quantitative description of membrane current and its application to conduction and excitation in nerve](#). *J Physiol* 117(4): 500-544.
9. Kiernan, M. C., Mogyoros, I., Hales, J. P., Gracies, J. M. and Burke, D. (1997). [Excitability changes in human cutaneous afferents induced by prolonged repetitive axonal activity](#). *J Physiol* 500 (Pt 1): 255-264.
10. Kress, G. J., Dowling, M. J., Meeks, J. P. and Mennerick, S. (2008). [High threshold, proximal initiation, and slow conduction velocity of action potentials in dentate granule neuron mossy fibers](#). *J Neurophysiol* 100(1): 281-291.
11. Kusano, K. (1966). [Electrical activity and structural correlates of giant nerve fibers in Kuruma shrimp \(*Penaeus japonicus*\)](#). *J Cell Physiol* 68(3): 361-383.
12. Meeks, J. P. and Mennerick, S. (2007). [Action potential initiation and propagation in CA3 pyramidal axons](#). *J Neurophysiol* 97(5): 3460-3472.
13. Miller, R. N. and Rinzel, J. (1981). [The dependence of impulse propagation speed on firing frequency, dispersion, for the Hodgkin-Huxley model](#). *Biophys J* 34(2): 227-259.
14. Rosenthal, J. J. and Bezanilla, F. (2000). [Seasonal variation in conduction velocity of action potentials in squid giant axons](#). *Biol Bull* 199(2): 135-143.
15. Smarandache, C. R. and Stein, W. (2007). [Sensory-induced modification of two motor patterns in the crab, *Cancer pagurus*](#). *J Exp Biol* 210(Pt 16): 2912-2922.
16. Städele, C. and Stein, W. (2016). [The site of spontaneous ectopic spike initiation facilitates signal integration in a sensory neuron](#). *J Neurosci* 36(25): 6718-6731.
17. Stein, R. B. and Pearson, K. G. (1971). [Predicted amplitude and form of action potentials recorded from unmyelinated nerve fibres](#). *J Theor Biol* 32(3): 539-558.
18. Swadlow, H. A. and Waxman, S. G. (1976). [Variations in conduction velocity and excitability following single and multiple impulses of visual callosal axons in the rabbit](#). *Exp Neurol* 53(1): 128-150.
19. Weidner, C., Schmelz, M., Schmidt, R., Hammarberg, B., Orstavik, K., Hilliges, M., Torebjork, H. E. and Handwerker, H. O. (2002). [Neural signal processing: the underestimated contribution of peripheral human C-fibers](#). *J Neurosci* 22(15): 6704-6712.