

Ovule Clearing Method for *Solanaceous* Species

Audrey Loubert-Hudon, Daniel P. Matton*

Institut de Recherche en Biologie Végétale, Département de Sciences Biologiques, Université de Montréal, Montréal, Québec, Canada

*For correspondence: dp.matton@umontreal.ca

[Abstract] Depending on the species, embryo sacs can be difficult to observe. Ovule clearing allows precise observation of whole tissues under Differential interference contrast (DIC) microscopy. The use of Methyl Salicylate as a clearing agent has proved to be particularly reliable for *Solanaceous* species.

Materials and Reagents

1. Flower buds
2. 37% Formaldehyde
3. Glacial acetic acid
4. Methyl Salicylate
5. Ethanol ($\geq 95\%$)
6. Fixative solution (FAA) (see Recipes)
7. Clearing solution (see Recipes)

Equipment

1. Forceps
2. Microscope with differential interference contrast (DIC)
3. Printed coated microscopy slides with wells (to avoid crushing the tissues)
4. Cover slips

Procedure

1. Carefully collect ovaries from flower buds or flowers at anthesis with precision forceps. During the process, keep ovaries and flowers on ice.
2. To facilitate penetration of the fixative and clearing solutions, the pericarp has to be removed. Under a magnifying glass, cut the top of the ovary's pericarp with precision forceps or scalpel to allow its removal (see Figure 2).
3. A bunch of ovules can then be carefully taken and fixed in FAA solution overnight. Gentle agitation at 4 °C helps to homogenize the fixation, but is not crucial.
4. The day after, fixed ovaries are incubated in 95-100% ethanol for at least one hour,

still with gentle agitation at room temperature.

5. The ovaries are then transferred in increasing ratio of Methyl Salicylate/EtOH solutions for 30 min each (1:3, 1:1, 3:1) at room temperature. Again, gentle agitation will help, but is not crucial.

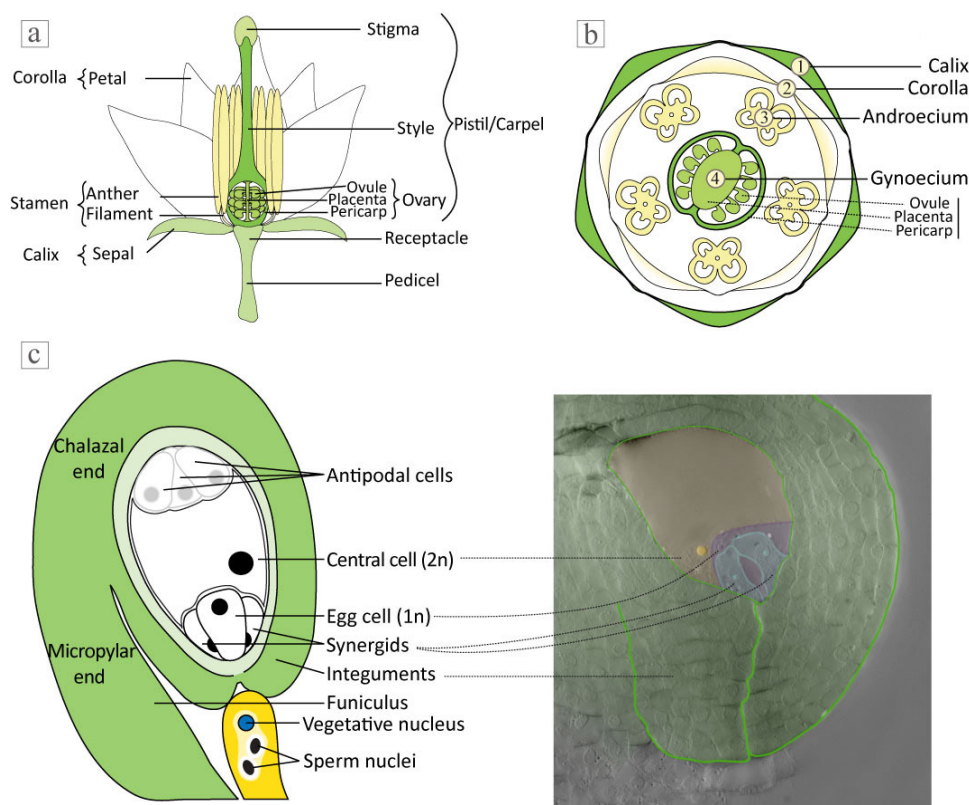


Figure 1. Organs and cells types of *Solanaceous* species' flower and female gametophyte. a. Longitudinal sections of a *Solanaceous* flower. b. Transversal sections of a *Solanaceous* flower. c. Description of ovule cells types and their visualization in cleared ovules (false colors are used to emphasize cells types).

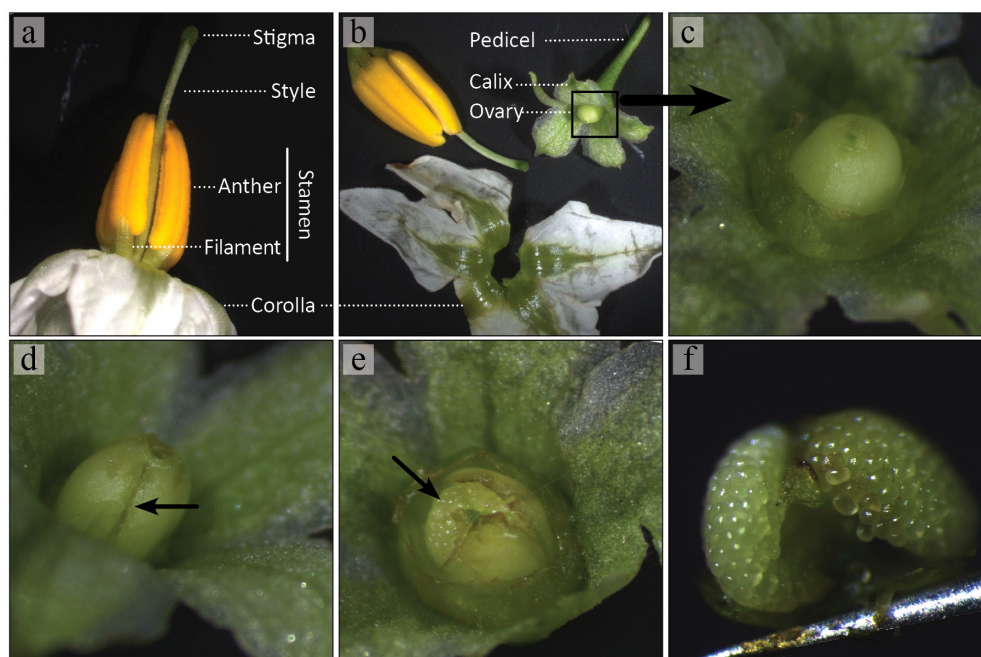


Figure 2. Steps by steps procedure to remove the pericarp of *Solanaceous* species. a. Collect flowers at the desired stage of development and keep them on ice; b. Remove the upper part of the flower carefully by hand to keep only the ovary still attached to the pedicel; c. d. Cut the top of the ovary's pericarp with precision forceps or scalpel to allow its removal; e. f. With the precision forceps, pinch under the bunch of ovules and put them directly in fixative solution.

6. The ovule can now be observed under differential interference contrast (DIC) microscopy. To avoid crushing the tissues, printed coated slides with wells are very useful. With precision forceps, carefully separate the ovules on the slide. Methyl Salicylate (100%) is used as mounting solution and slides are sealed with clear nail polish.

Recipes

1. Fixative solution (FAA)

1% formaldehyde

0.5% glacial acetic acid

50% ethanol

Note: It can be prepared in advance and stored at room temperature.

2. Clearing solution (see Procedure step 5, for concentration/ratio)

Methyl salicylate

Ethanol

Acknowledgments

This protocol is adapted from Chevalier *et al.* (2013) and Estrada-Luna *et al.* (2004).

References

1. Chevalier, E., Loubert-Hudon, A. and Matton, D. P. (2013). [ScRALF3, a secreted RALF-like peptide involved in cell-cell communication between the sporophyte and the female gametophyte in a *solanaceous* species](#). *Plant J* 73(6): 1019-1033.
2. Estrada-Luna, A., Garcia-Aguilar, M. and Vielle-Calzada, J. P. (2004). [Female reproductive development and pollen tube growth in diploid genotypes of *Solanum cardiophyllum* Lindl.](#) *Sexual Plant Reproduction* 17(3): 117-124.